



LOFMEM™ T PTFE membrane filter cartridges

Eaton's LOFMEM T membrane filter cartridges are suitable for gas/vent applications and the filtration of aggressive compounds. Specific uses include: strong acids/bases, compressed gases, photoresists, pharmaceutical intermediates, hot DI water and fermentation air.

LOFMEM T membrane filter cartridges offer precision filtration – even under harsh, corrosive environments. The polytetrafluoroethylene (PTFE) membrane with additional polypropylene support layers and components offer superior hydrophobicity and water intrusion resistance compared to PVDF-based media and are more economical than all fluorocarbon filters.

Features and benefits

- High flow rates and high surface area minimizes total system size requirements
- Full traceability marking
- Integrity testing and 100% flushing with 18 MΩ-cm de-ionized water prior to shipment
- Manufactured in ISO Class 7 Cleanroom Environment

Design

Filter material
PTFE membrane

Inner core, cage, end caps
Polypropylene

Support layers
Polypropylene

O-rings
Silicone (standard), Buna-N, EPDM, FPM, FEP encapsulated FPM

Retention ratings
0.05, 0.1, 0.2, 0.45, 1 µm

Technical data

Nominal lengths
10", 20", 30", 40"
(254, 508, 762, 1016 mm)

Outside diameter
2.7" (69 mm)

Inside diameter
1" (25.4 mm)

Surface area
8.5 ft² (0.79 m²) per 10" element

Max. operating temperature
203 °F (95 °C)

Max. differential pressures forward
75 psid @ 70 °F (5.2 bar @ 21 °C)
40 psid @ 176 °F (2.8 bar @ 80 °C)
15 psid @ 203 °F (1.0 bar @ 95 °C)

Max. differential pressure reverse
40 psid @ 70 °F (2.8 bar @ 21 °C)




Powering Business Worldwide

LOFMEM T membrane filter cartridges

Flow rate*

(68 °F/20 °C per 10" filter cartridge for water)

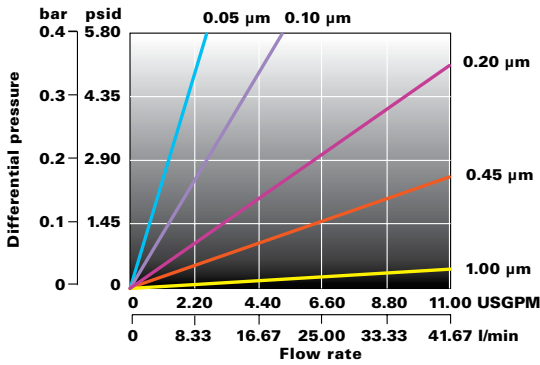

* For liquids other than water, multiply pressure drop by fluid viscosity in centipoise.

Performance specifications

Steam/autoclave


Filter cartridges will withstand at least 50 steam/autoclave 1 hour cycles @ 135 °C

Integrity test values per 10" filter cartridge

Air Diffusion per 10-inch filter cartridge wet with 60/40 IPA/water. Contact Eaton for specific method.

Pore size	Specification
0.05 µm:	≤ 50 cc/min @ 1.5 bar
0.1 µm:	≤ 50 cc/min @ 1.2 bar
0.2 µm:	≤ 20 cc/min @ 0.8 bar
0.45 µm:	≤ 15 cc/min @ 0.34 bar
1 µm:	≤ 15 cc/min @ 0.2 bar

Ordering code

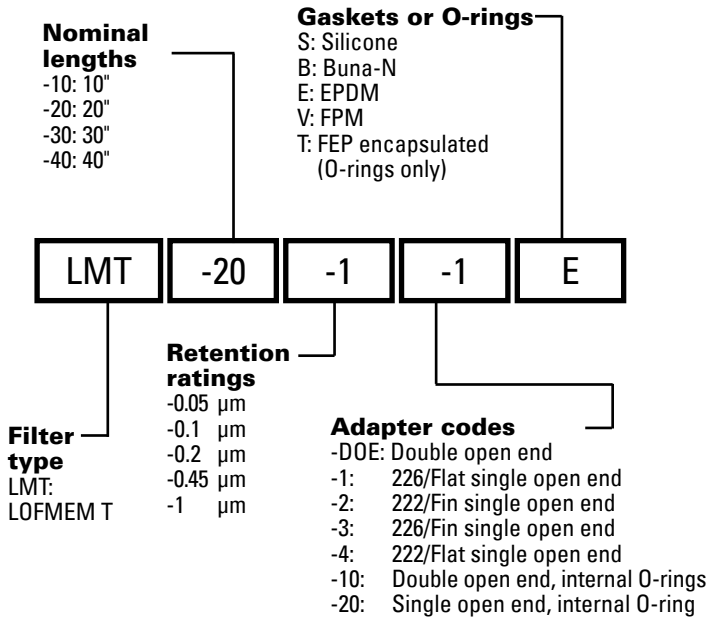

North America
 44 Apple Street
 Tinton Falls, NJ 07724
 Toll Free: 800 656-3344
 (North America only)
 Tel: +1 732 212-4700

China
 No. 3, Lane 280,
 Linhong Road
 Changning District, 200335
 Shanghai, P.R. China
 Tel: +86 21 5200-0099

Europe/Africa/Middle East
 Auf der Heide 2
 53947 Nettersheim, Germany
 Tel: +49 2486 809-0

Singapore
 100G Pasir Panjang Road #07-08
 Singapore 118523
 Tel: +65 6825-1668

Friedensstraße 41
 68804 Altlußheim, Germany
 Tel: +49 6205 2094-0
 An den Nahewiesen 24
 55450 Langenlonsheim, Germany
 Tel: +49 6704 204-0

Brazil
 Rua Clark, 2061 - Macuco
 13279-400 - Valinhos, Brazil
 Tel: +55 11 3616-8400

For more information, please email us at filtration@eaton.com or visit www.eaton.com/filtration

© 2018 Eaton. All rights reserved. All trademarks and registered trademarks are the property of their respective owners. All information and recommendations appearing in this brochure concerning the use of products described herein are based on tests believed to be reliable. However, it is the user's responsibility to determine the suitability for his own use of such products. Since the actual use by others is beyond our control, no guarantee, expressed or implied, is made by Eaton as to the effects of such use or the results to be obtained. Eaton assumes no liability arising out of the use by others of such products. Nor is the information herein to be construed as absolutely complete, since additional information may be necessary or desirable when particular or exceptional conditions or circumstances exist or because of applicable laws or government regulations.

EN
 EF-LMT
 05-2018