Enclosed Plate and Frame Filter BECO INTEGRA® PLATE 600 EP

Stainless Steel Chassis, Plastic Filter Plates, and Cake Frames

BECO INTEGRA PLATE 600 EP enclosed plate and frame filter is the depth filtration system for critical filtration tasks in the pharmaceutical and biotech, chemical, cosmetics, and food industries.

BECO INTEGRA PLATE 600 EP is a multi-sheet filter with external chamber. Depth filter sheets seal the product channels. Gaskets at the external chamber seal the entire filter package.

Special Features

- CIP of filter package without depth filter sheets.
- **Complete drainage** of filter due to unique filter plate design.
- **Safe operation:** No drip loss due to sealed external chamber and automatic hydraulic pressure regulation.
- Easy handling and good cleanability thanks to ergonomic design of chassis, and open support grid of filter plates; cake discharge downwards into tray underneath of filter package.
- **High flexibility:** The filter can be used as a sheet filter or as a pre-coat filter. For pre-coating/removing of high particle load, cake frames in different widths are available.
- BECO INTEGRA PLATE 600 EP filter is designed and manufactured according to the cGMP guidelines and the FDA regulations.

Customization or additional components allow for system expansion or complete systems according to customer specifications.

Technical Data

Operating pressure/ temperature:	max. 87 psi (600 kPa/6 bar) at 104 °F (40 °C) max. 43.5 psi (300 kPa/3 bar) at 140 °F (60 °C)
Differential pressure/ temperature:	max. 43.5 psi (300 kPa/3 bar) at 104 °F (40 °C)
Filter area per sheet:	2.9 ft² (0.27 m²)





Materials

Product contacting parts:	Polypropylene/ Stainless steel AISI 316L
Other parts:	Stainless steel AISI 304
Gaskets/O-rings:	Silicone, EPDM (FDA listed), Viton



Туре	Filter a	area	Cak volun		Chassis size	Filter plates/ cake frames	Depth filter sheets 20.1 x 24.8 in (512 x 630.5 mm)	We	eight
	ft²	(m²)	gal	(I)		qty.	qty.	≈ Ib	(kg)
600 SF	75.6	(7.02)	14.0	(53)	1	12/13	26	2,976	(1,350)
with 0.79"	133.7	(12.42)	24.6	(93)	2	21/22	46	3,307	(1,500)
(20 mm)	191.8	(17.82)	35.4	(134)	3	32/33	66	3,693	(1,675)
cake frames	32.3	(27.00)	53.4	(202)	4	49/50	100	4,321	(1,960)
600 ASF	64.0	(5.94)	19.6	(74)	1	10/11	22	2,921	(1,325)
with 1.18"	110.4	(10.26)	33.1	(128)	2	18/19	38	3,285	(1,490)
(30 mm)	162.75	(15.12)	50.0	(189)	3	27/28	56	3,660	(1,660)
cake frames	244.1	(22.68)	74.8	(283)	4	41/42	84	4,277	(1,940)
600 ASF	52.3	(4.86)	22.5	(85)	1	8/9	18	2,943	(1,335)
with 1.57"	98.8	(9.18)	42.5	(161)	2	16/17	34	3,307	(1,500)
(40 mm)	139.5	(12.96)	60.0	(227)	3	23/24	48	3,638	(1,650)
cake frames	209.3	(19.44)	89.8	(340)	4	31/32	72	4,101	(1,860)
600 ASF	40.7	(3.78)	29.6	(103)	1	6/7	14	2,932	(1,330)
with 2.36"	75.6	(7.02)	55.0	(193)	2	12/13	26	3,263	(1,480)
(60 mm)	110.4	(10.26)	80.3	(282)	3	18/19	38	3,594	(1,630)
cake frames		(15.12)	118.4	(415)	4	27/28	56	4,145	(1,880)

* usable cake volume

Filter Plates/Cake Frames


Width	Cake volume*	Weight	
in (mm)	gal (I)	lb (kg)	
1.02 (26)		19.9 (9.0)	
1.02 (26)		17.6 (8.0)	
0.79 (20)	1.07 (4.05)	6.6 (3.0)	
1.18 (30)	1.78 (6.75)	9.9 (4.5)	
1.57 (40)	2.50 (9.45)	13.2 (6.0)	
2.36 (60)	3.92 (14.85)	19.9 (9.0)	
	in (mm) 1.02 (26) 1.02 (26) 0.79 (20) 1.18 (30) 1.57 (40)	volume* in (mm) gal (l) 1.02 (26) - - 0.79 (20) 1.07 (4.05) 1.18 (30) 1.78 (6.75) 1.57 (40) 2.50 (9.45)	

Connections

Filter inlet:	TC 11/2" ISO2852/ASME BSP
Filter outlet:	TC 11/2" ISO2852/ASME BSP
External chamber:	TC 11/2" ISO2852/ASME BSP

Other Connections Threaded connections according to DIN 1151/11864

1445

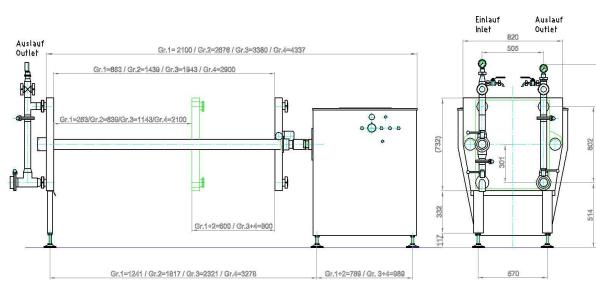

usable cake volume

North America 44 Apple Street Tinton Falls, NJ 07724 Toll Free: 800 656-3344 (North America only) Tel: +1 732 212-4700

Europe/Africa/Middle East

Auf der Heide 2 53947 Nettersheim, Germany Tel: +49 2486 809-0

Friedensstraße 41 68804 Altlußheim, Germany Tel: +49 6205 2094-0

An den Nahewiesen 24 55450 Langenlonsheim, Germany Tel: +49 6704 204-0

China

No. 3, Lane 280, Linhong Road Changning District, 200335 Shanghai, P.R. China Tel: +86 21 5200-0099

Singapore

100G Pasir Panjang Road #07-08 Singapore 118523 Tel: +65 6825-1668

Brazil

Av. Ermano Marchetti, 1435 -Água Branca, São Paulo - SP, 05038-001, Brazil Tel: +55 11 3616-8461

For more information, please email us at filtration@eaton.com or visit www.eaton.com/filtration

ΕN 1 A 2.6.11.6 08-2019

© 2019 Eaton. All rights reserved. All trademarks and registered trademarks are the property of their respective owners. All information and recommendations appearing in this brochure concerning the use of products described herein are based on tests believed to be reliable. However, it is the user's responsibility to determine the suitability for his own use of such products. Since the actual use by others is beyond our control, no guarantee, expressed or implied, is made by Eaton as to the effects of such use or the results to be obtained. Eaton assumes no liability arising out of the use by others of such products. Nor is the information herein to be construed as absolutely complete,

since additional information may be necessary or desirable when particular or exceptional conditions or circumstances exist or because of applicable laws or government regulations.

